


## Instant Notifier Plus Workstation Web App

Instant Notifier Plus is a web app that users can access through the Freezer Monitoring web application to receive audio alerts on their computers. This web app is helpful for receiving real-time audio alerts while using your workstation during work hours.

There are two ways to access the Instant Notifier Plus workstation web app: through the Freezer Monitoring web app or by adding the Instant Notifier Plus extension to your browser.


### A. Through the Freezer Monitoring Web App:

1. Log in to the freezer monitoring application. On the welcome page, you will see “My Applications,” including:
  - **Environmental Monitoring** (to check the status of freezer monitors)
  - **Instant Notifier Plus** (to receive alert notifications)

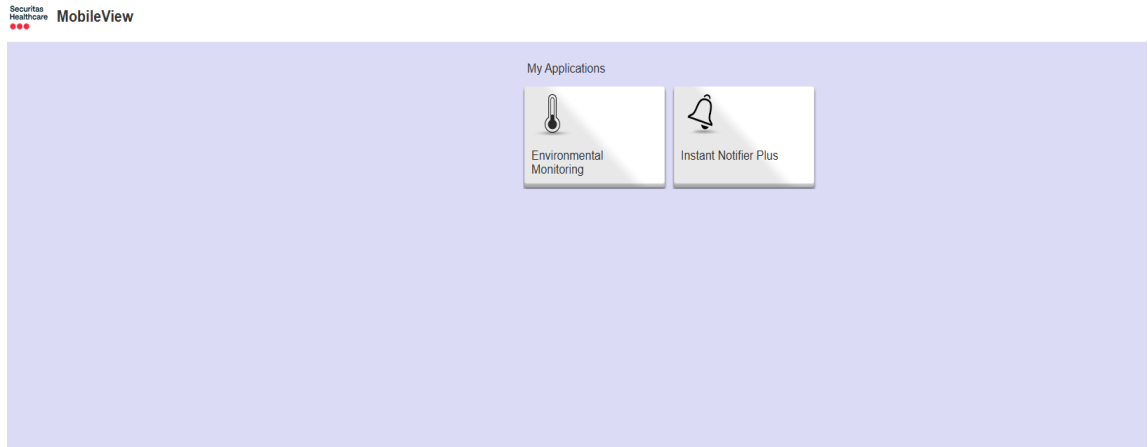


**Note:** You can open Environmental Monitoring, Instant Notifier Plus, or both. To open both, right-click on one application and select “Open link in a new tab,” then open the other in a separate tab.


Another way to switch between applications is by using the drop-down arrow in the top-left corner after selecting an application. You can also right-click and select “Open link in a new tab” to open another application. This also applies if you have administrative access and want to switch between applications on the administration page.

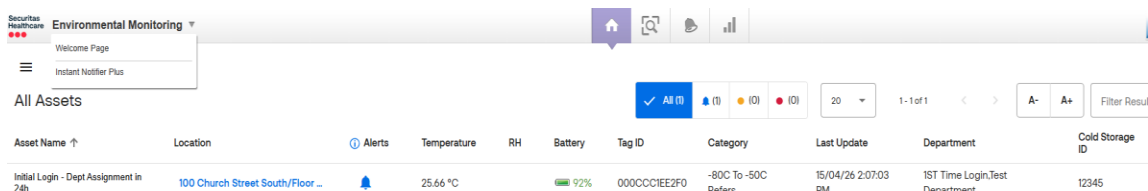


- After selecting the Instant Notifier Plus application, toggle the speaker icon on (top-left corner) to enable audio alerts. To turn off audio alerts, toggle the speaker icon off and you can always close the tab to stop the alerts.

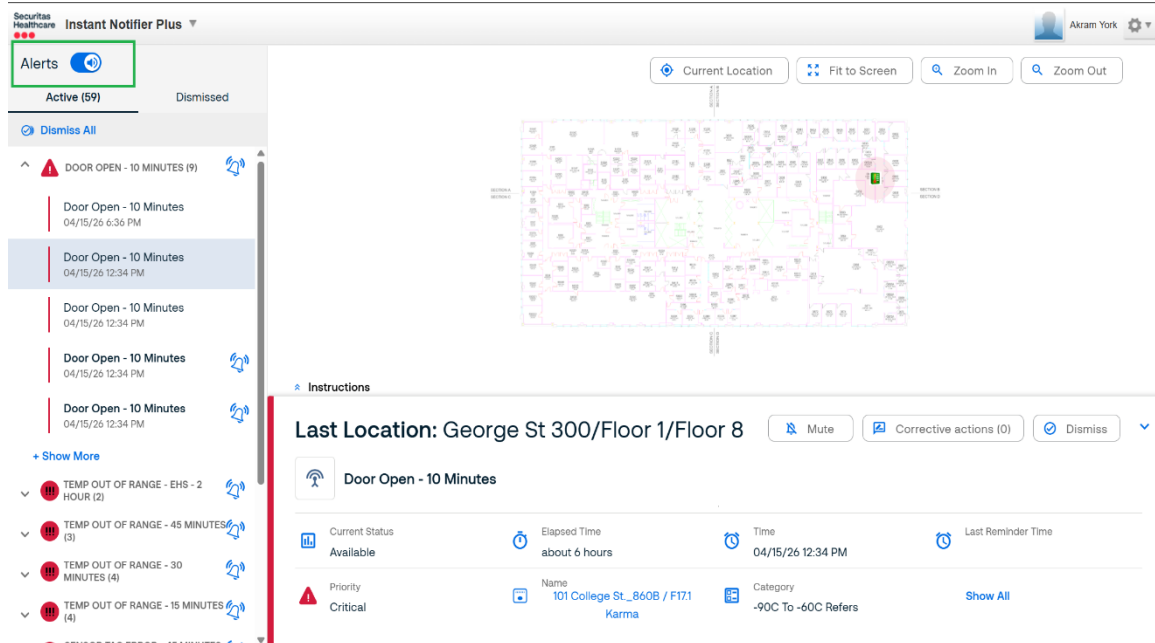

- To mute a single alert, click the **Mute** button in the middle of the page.

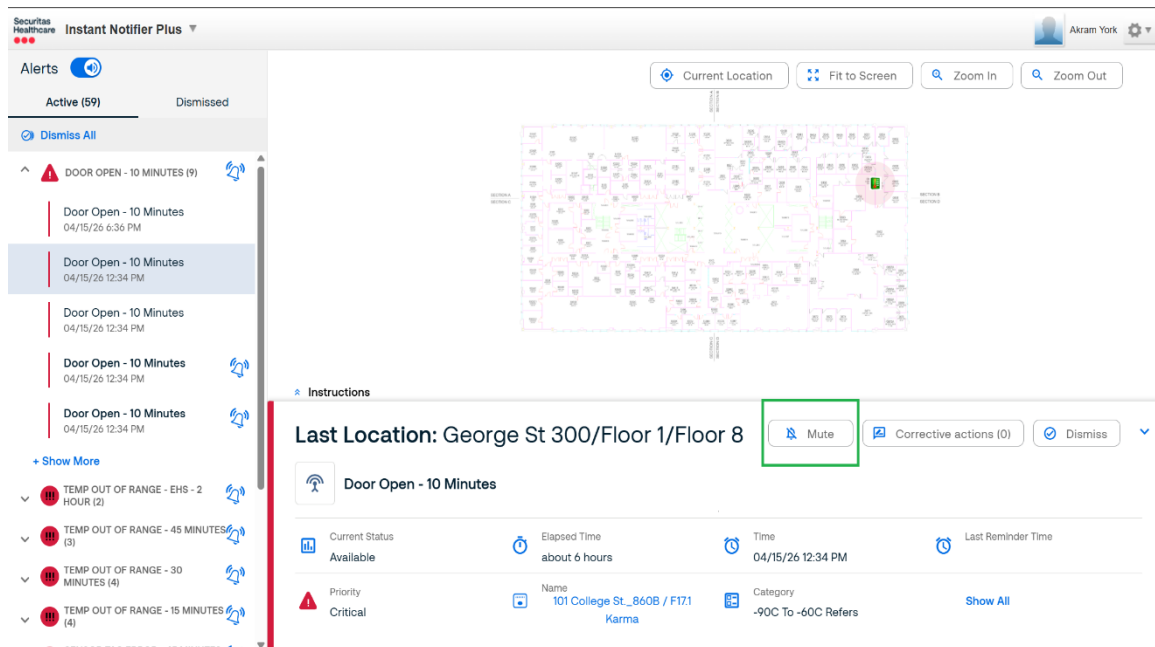


- To dismiss (acknowledge) an alert, click the **Dismiss** button. A new window will appear prompting you to add Corrective Actions, similar to dismissing alerts in the Environmental Monitoring application.

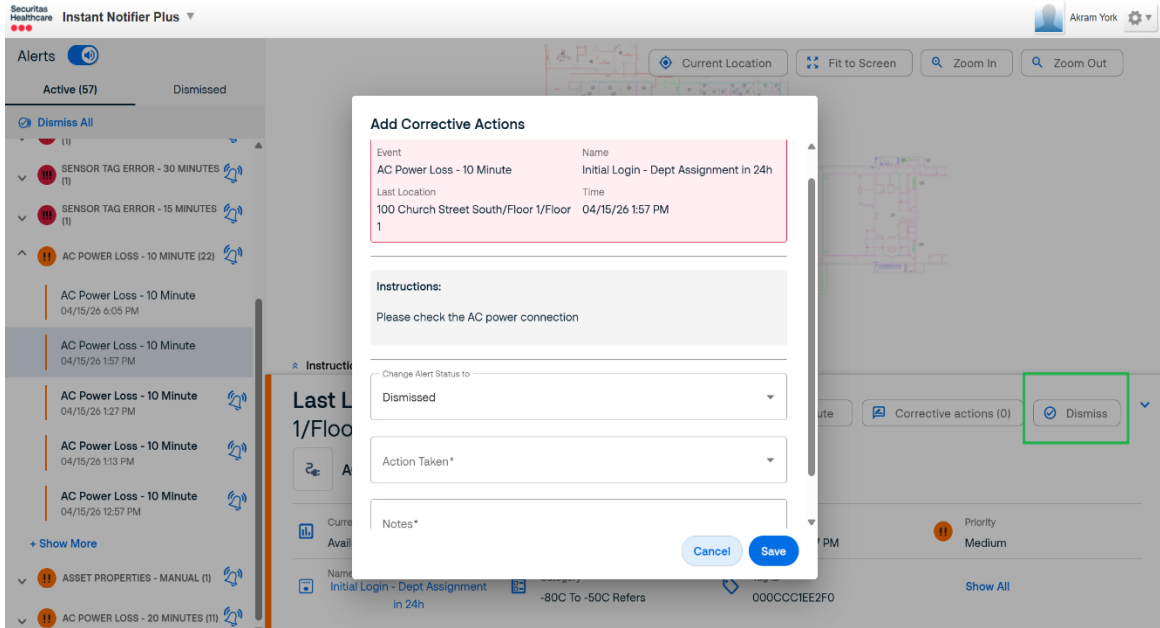


## B. Through the Instant Notifier Plus Browser Extension:


- Search for “**Instant Notifier Plus**” in your browser’s extension store and add it to your browser. For example, Chrome users can install it from the [Chrome Web Store](#) link directly.

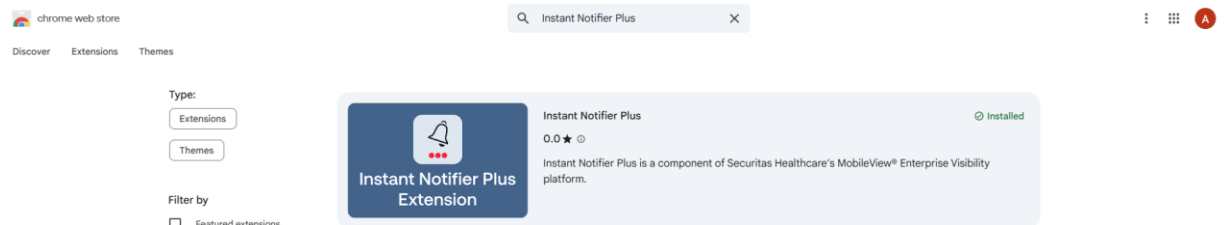


- After installation, open the extension and configure it by entering the server name: <https://freezermonitor.its.yale.edu>, along with your NetID and password (SSO login is not supported in the extension at this time).

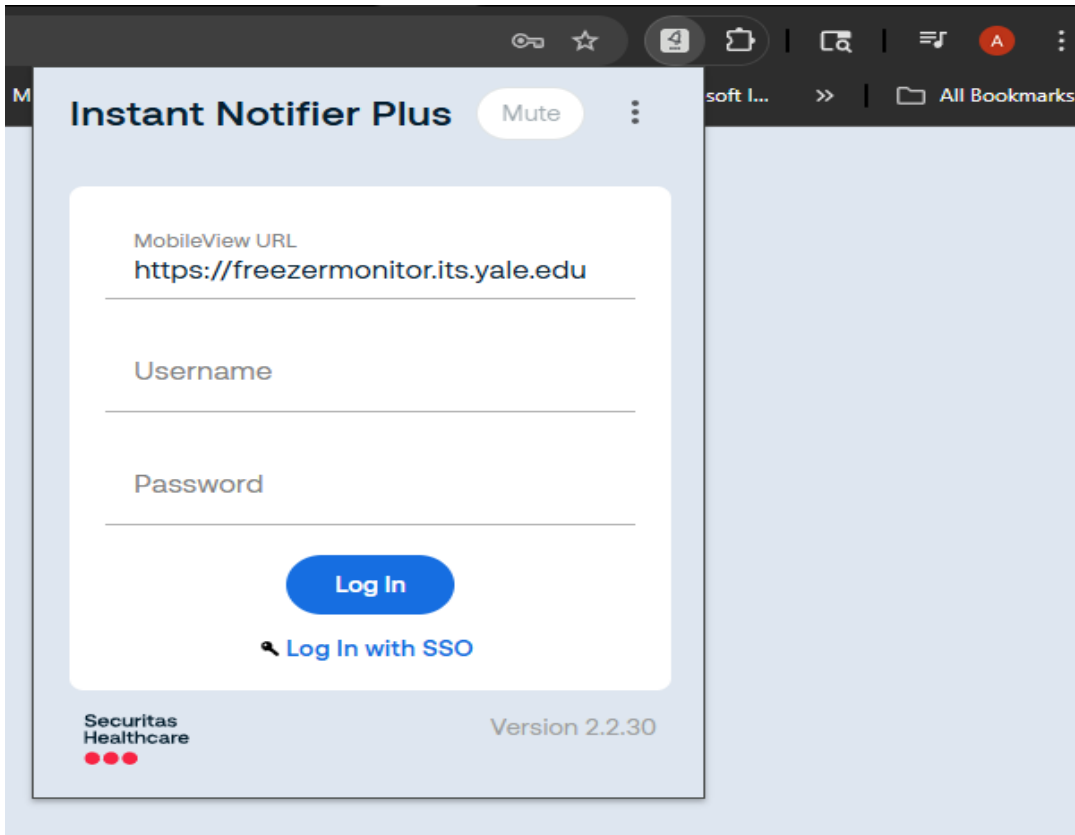

3. Once logged in, you will begin receiving recent and new alerts. The Instant Notifier Plus application page will also pop up to display alert details.

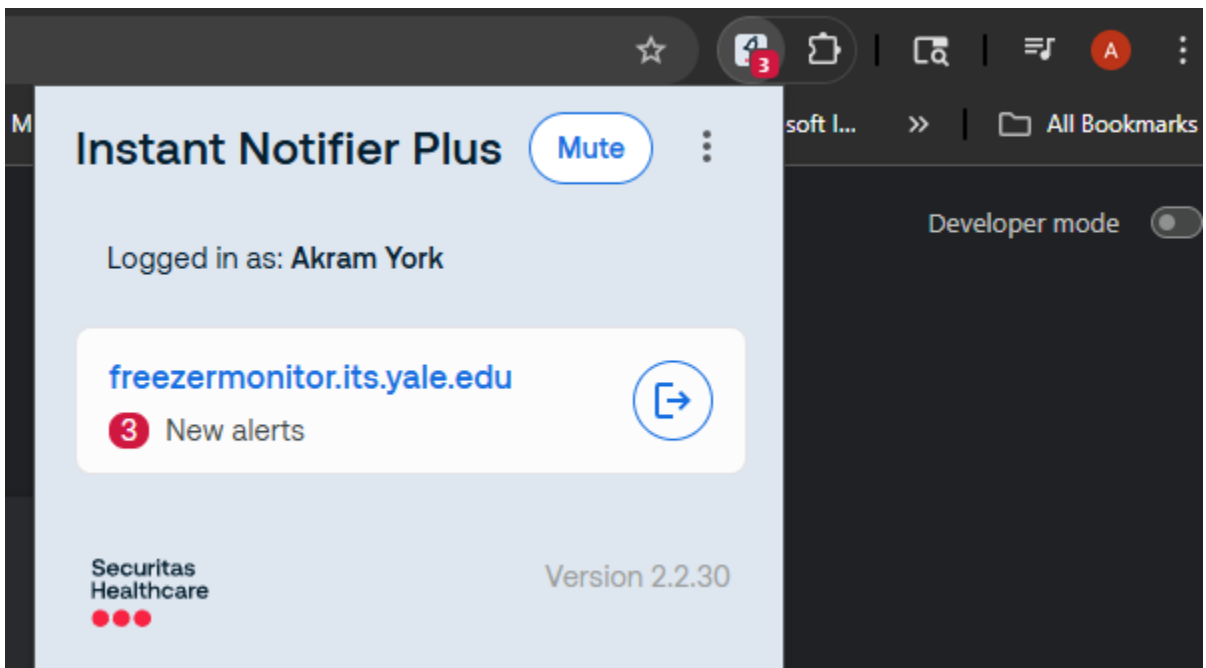

- When you receive an alert, you can mute it by clicking the **Mute** button. To dismiss the alert and add a corrective action, use the pop-up Instant Notifier Plus application (as described above).

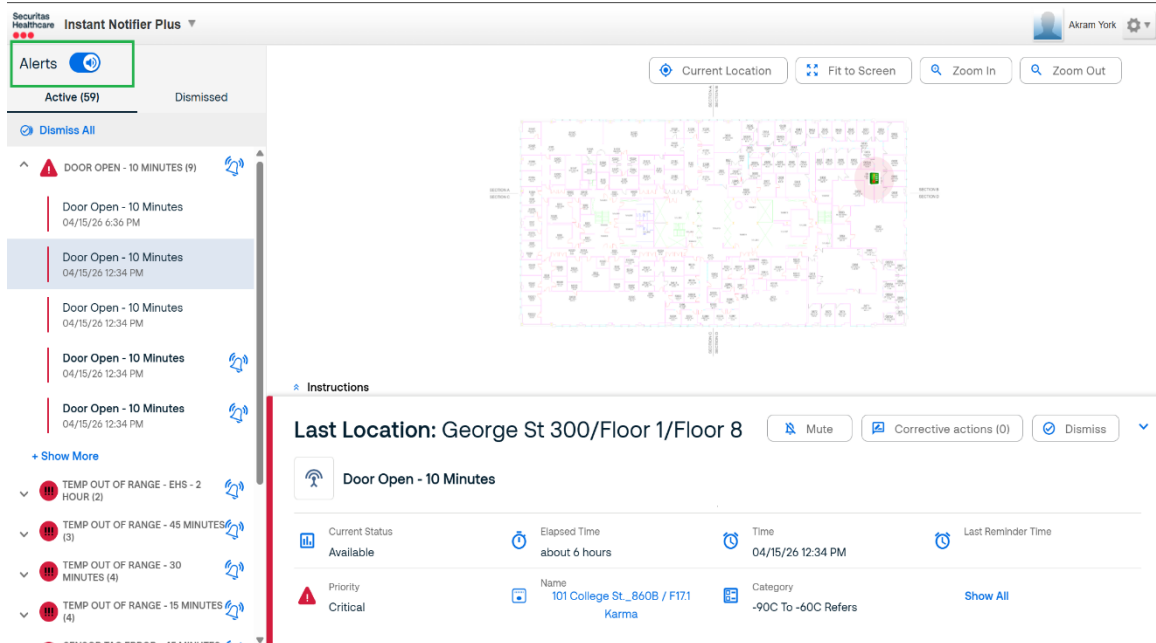

- You can log out of the extension at any time by clicking the Exit icon.

